




UNITE[®]
FOOT & ANKLE

Soft Tissue Fixation Solutions
Intelligently designed for tendon and ligament repair.


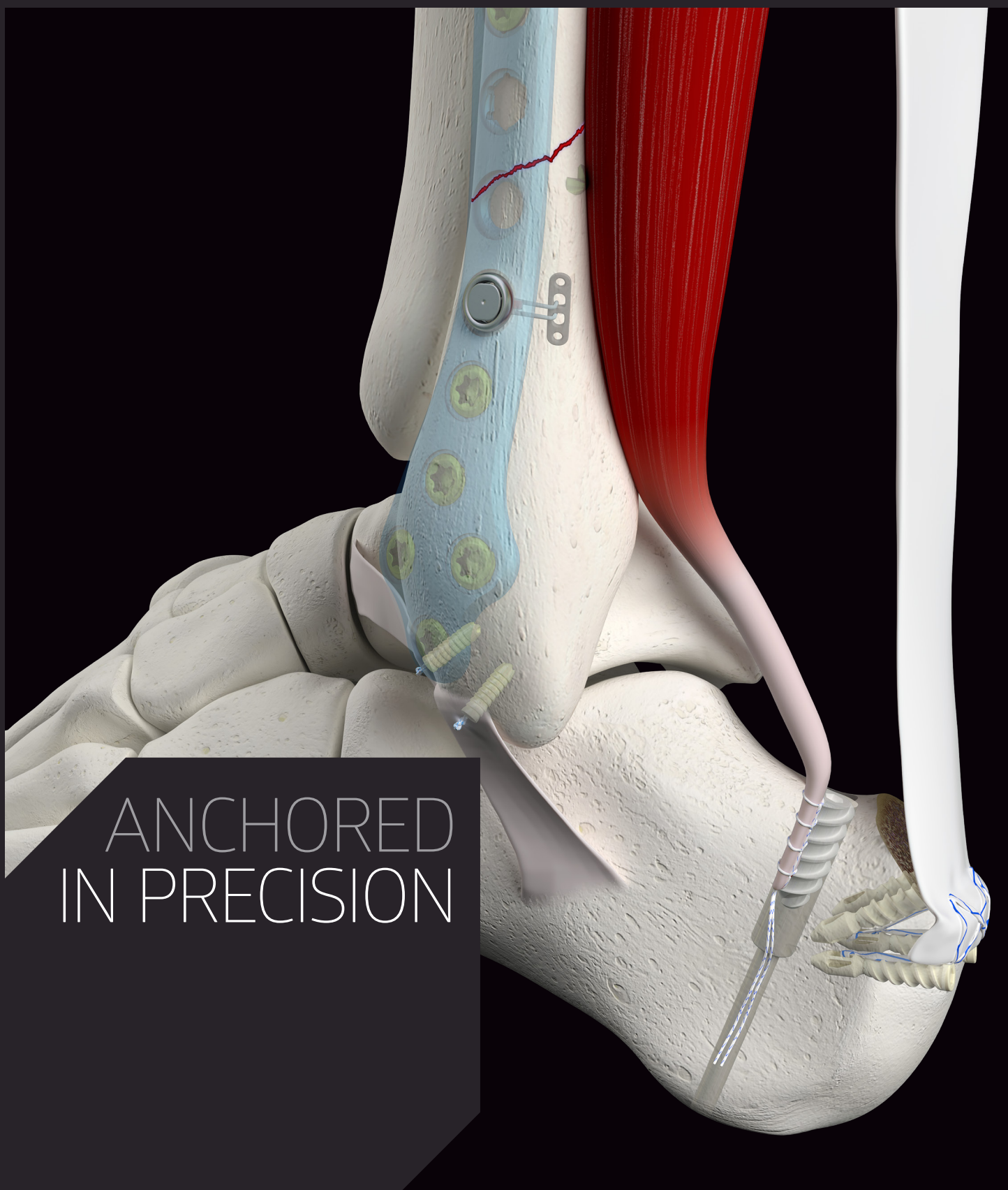


ANCHORED
IN PRECISION

SYNDEX™

WITH
CONSTRUCTOR®
TECHNOLOGY

- A knotless adjustable button technology indicated for syndesmosis repair
- Features a novel, self-locking design to help prevent loss of reduction
- Superior fixation under cyclic loading
- Substantial difference in load-to-failure compared to the market leader¹
- Button fits within UNITE Lateral Fibula Plate with minimal prominence

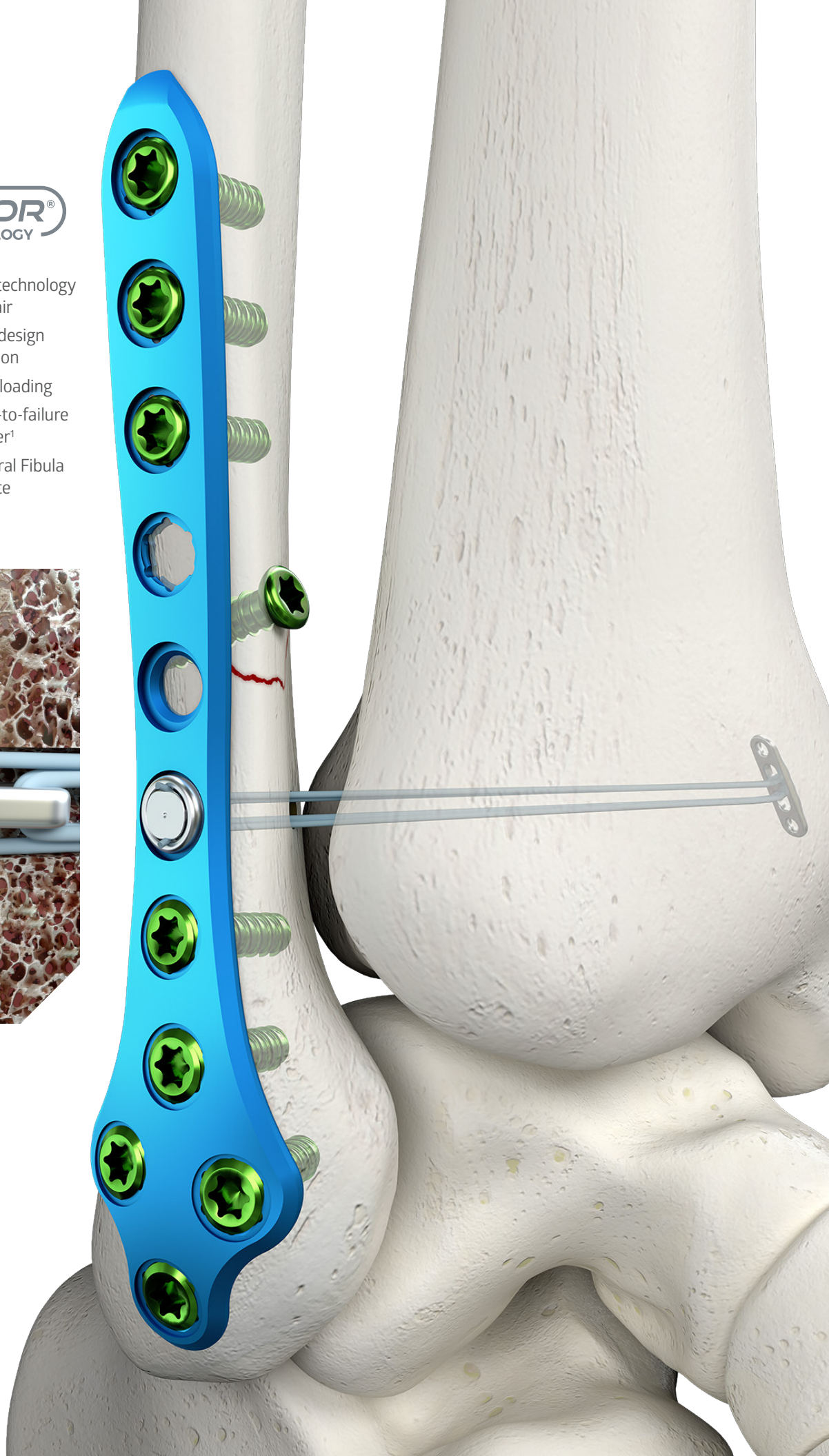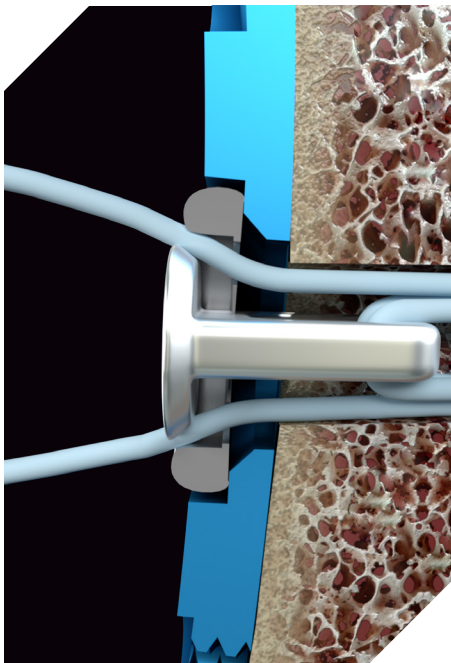


DEXLOCK™ KNOTLESS Suture Anchors

- Made from PEEK with an anchor design that allows for suture or tape tensioning to the desired repair location after the eyelet is fully seated
- Two-hole eyelet design accommodates multiple sutures or tapes for double row knotless repairs
- Dual-lead thread design promotes faster insertion
- Anchors provide 500 N of pull-out strength
- Ultra-high molecular weight polyethylene tape is less abrasive than a leading competitor¹
- Tape features round-to-flat braid construction for easy handling and passing; lies flat with no inner core to distribute a consistent, even load over a greater surface area

Sizes available:

- Ø4.5 x 20mm
- Ø5.5 x 20mm

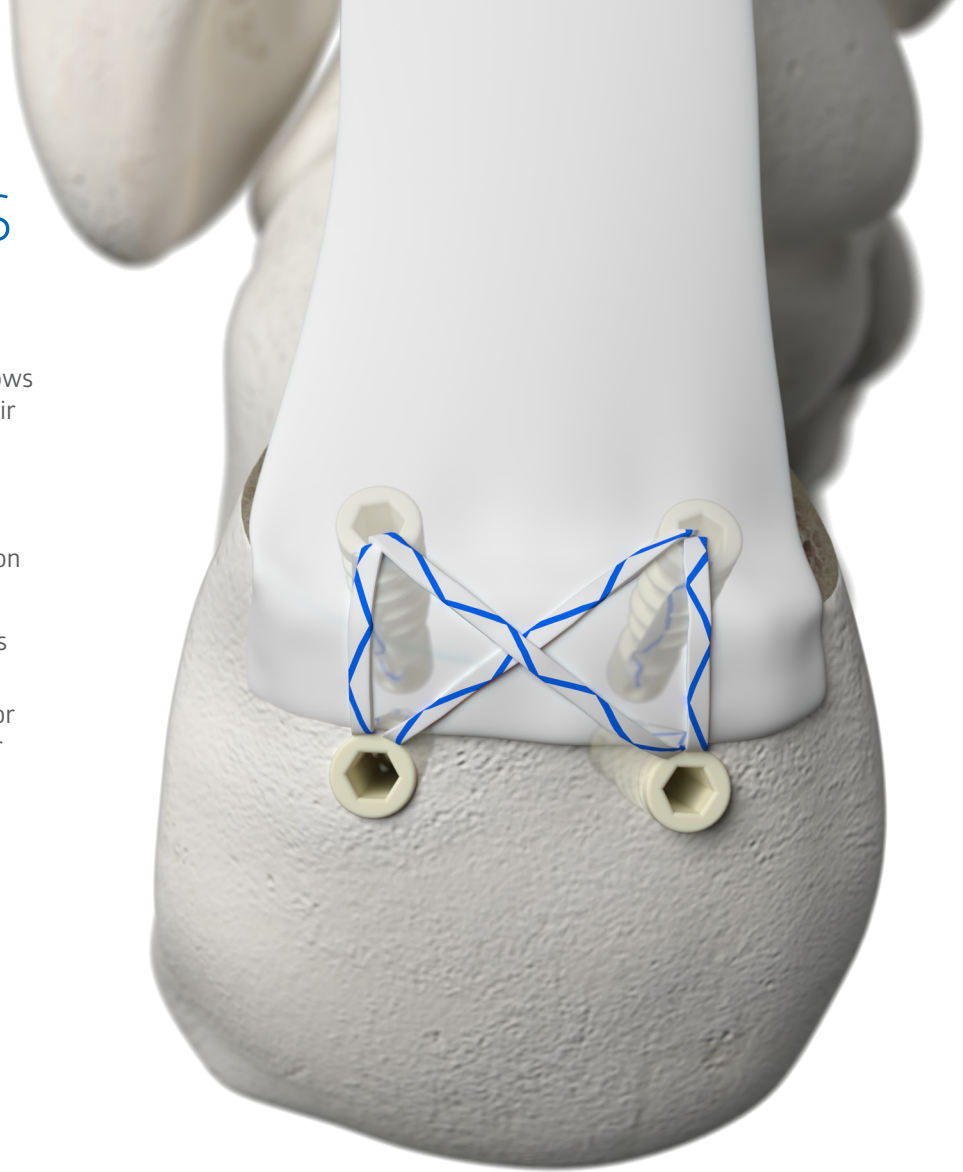


DEXTACK™ PUSH Suture Anchors

- Made from PEEK and double loaded with non-absorbable, ultra-high molecular weight polyethylene suture
- Compact anchor body design is ideal for lateral ankle ligament (ATFL) repair, medial deltoid ligament repair, and kidner procedures
- Approximately 300 N in pull-out strength

Sizes available:

- Ø2.9 x 10mm, (2) #0 suture, (4) CT-2 needles
- Ø3.3 x 10mm, (2) #2 suture, (4) MO-6 needles

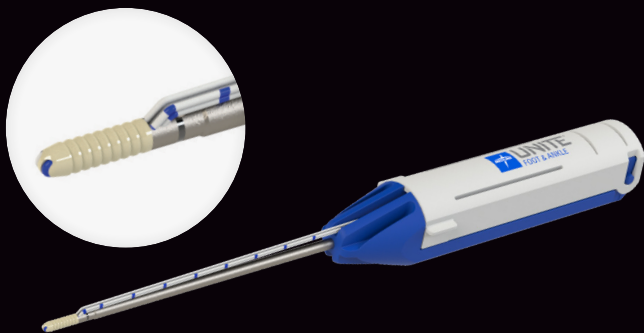



DEXTACK TWIST Suture Anchors

- Made from PEEK and double loaded with non-absorbable, ultra-high molecular weight polyethylene suture
- Dual-lead thread pitch design is ideal for achilles repair procedures with approximately 300 N of pull-out strength

Sizes available:

- Ø4.5 x 15mm, (2) #2 suture, (4) MO-6 needles
- Ø5.5 x 15mm, (2) #2 suture, (4) MO-6 needles

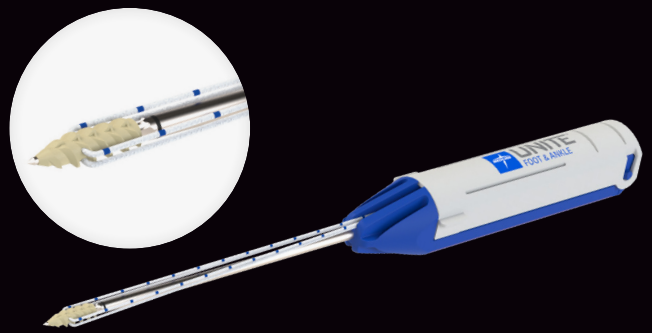

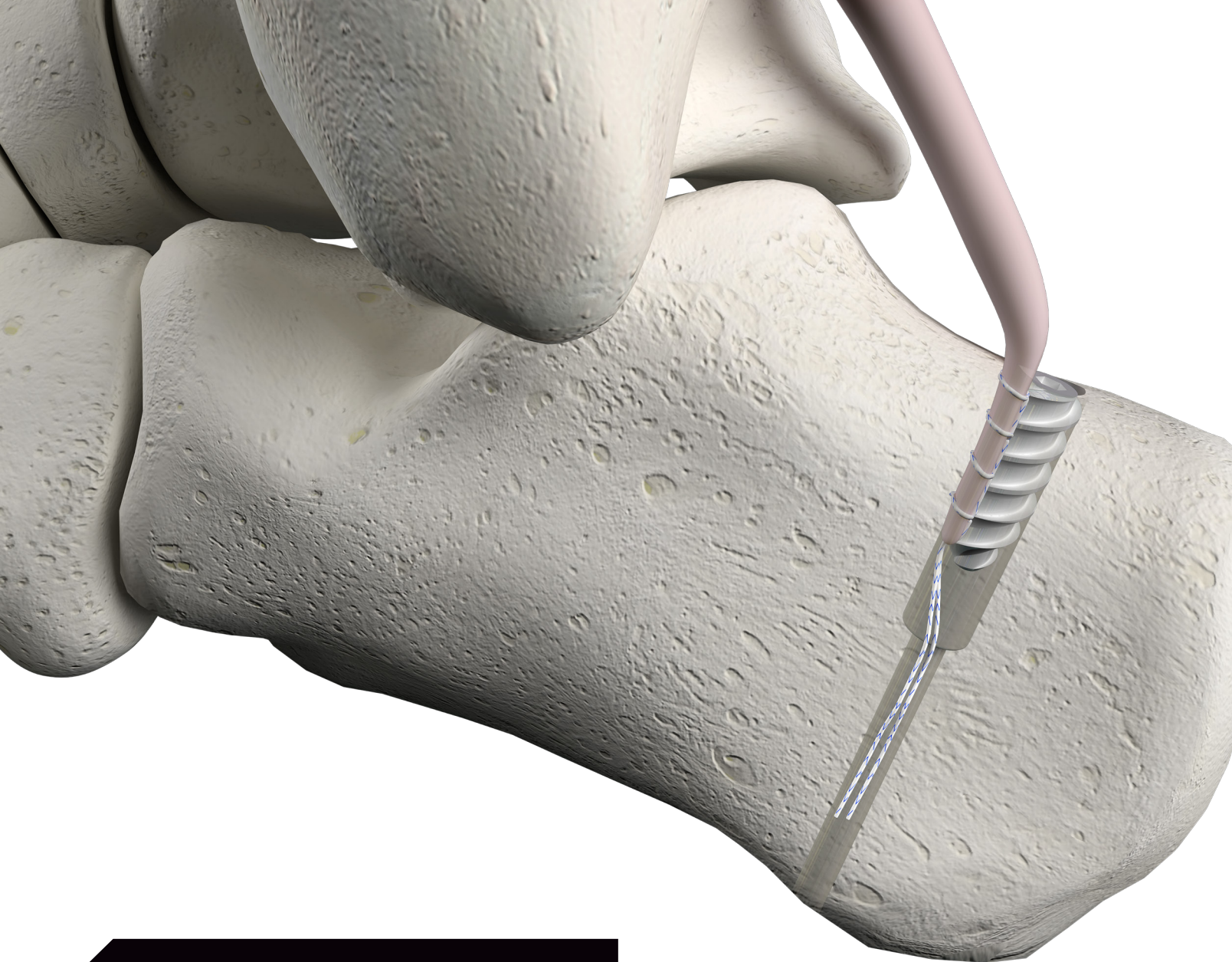


DEXTEN™ Tenodesis Screws

- Made from PEEK
- Available implant diameters and lengths are ideal for common tendon transfers in foot and ankle surgery, such as FDL and FHL tendon transfers
- Includes suture loop to whipstitch tendons, a tendon sizer, **Ø2.4 x 240mm suture guide drill pins**, acorn reamers, and a hex driver

Sizes available:

- Ø4.5 x 15mm
- Ø5.0 x 15mm
- Ø5.5 x 15mm
- Ø6.0 x 15mm
- Ø6.5 x 15mm
- Ø7.0 x 15mm

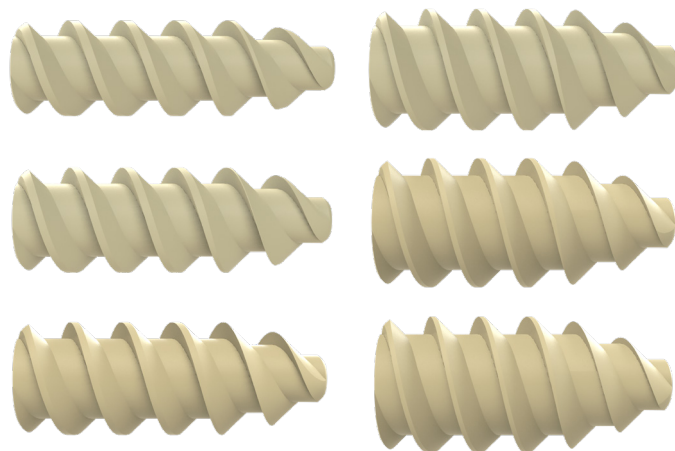

Optimized for Foot & Ankle

Properly sized implant options, shorter instrumentation, and convenient kit configurations. The Medline UNITE Soft Tissue Fixation portfolio executes on our philosophy of intelligent design with products tailored to address the unique needs of the Foot & Ankle specialist.



To learn more or schedule a case, contact your Medline UNITE Representative or visit medlineunite.com.



Medline Industries, LP Three Lakes Drive, Northfield, IL 60093
Medline United States
1-800-MEDLINE (633-5463)
medlineunite.com | unite@medline.com

Medline Canada
1-800-268-2848 | medline.ca | canada@medline.com
Medline México
01-800-831-0898 | medlinemexico.com | mexico@medline.com

Follow us [f](#) [in](#) [t](#) [BLOG](#)

REFERENCES. 1. Data on file

We reserve the right to correct any errors that may occur within this brochure.

© 2022 Medline Industries, LP. All rights reserved. Constrictor® is a registered trademark of Dunamis Medical. Medline is a registered trademark of Medline Industries, LP.
MKT19W2320989 / LIT997 / 1M / SG / 78